



OXYGEN UNLIMITED

NUVO

Oxygen concentrator

High-efficiency

Cost-effective

Reliable



Nidek Medical's Nuvo Oxygen Concentrator carries forward our commitment of providing reliable, cost-effective and quality products.


An Oxygen Concentrator is a machine that generates medical oxygen. It runs on electricity and the only maintenance required is the regular cleaning of filters.

The Nuvo Oxygen Concentrator has been designed to provide a combination of enhanced features, reliability of a proven technology and unsurpassed quality.

The Nuvo Oxygen Concentrator is the largest selling concentrator in India. We promise you complete satisfaction and convenience.

- Quiet Operation - less than 48 dBA
- Lightweight - 50 lbs. - Easy to Transport and Maneuver
- Slender & Sleek Cabinet Design with Integrated Handles
- Quick Snap Rear Panel allows easy access to filter, hourmeter & battery
- Superior Grade of Molecular Sieve
- Patented Maintenance Free Rotary Poppet Valve
- Exclusive Moisture Reducing Wick Assembly
- Reinforced Base Substructure and Recessed Casters to guard against breakage
- Spring Mounted High Efficiency Thomas Nexus Compressor
- Patented RPSA Technology

Largest selling concentrator in India

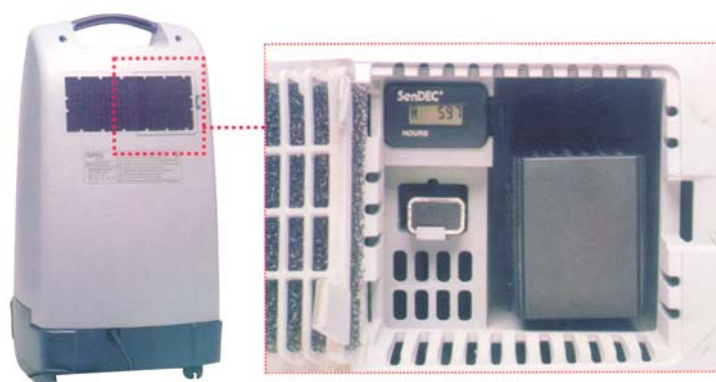

Specifications

Size (LxWxH)	394 X 396 X 706 mm
Weight	22.7 Kgs
Sound Level	< 47 dBA
Power Consumption	400 Watts average
Safety Alarms	Power Failure, Process Failure, High & Low Pressure Low Oxygen Purity (OMS model)
Oxygen Purity	0-2 LPM - 94%, 2-4 LPM - 93%, 5 LPM - 90%
Output Pressure	7 psi
Oxygen Output	Variable - 0 to 5 LPM (1/2 Liter increment)
Filters	Cabinet, Compressor Intake & Bacteria
Warranty	One Year



Nidek Medical India Private Limited

Head Office: 238B A.J.C. Bose Road, Kolkata 700020
Tel: (033) 4017 7000 (30 Lines), Fax: 91 (033) 2289 4217

Delhi: 238 Ansal Chambers II, 6, Bhikaji Cama Palace, New Delhi 110066
Tel: (011) 2616 5742, 93501 83222, 93135 774444

Chennai: AE - 32, (Old No 17), 7th Street, 10th Main Road, Anna Nagar, Chennai 600040
Tel: (044) 4266 3662, 93808 06777

Mumbai: 211 B, Damji-Shamji Industrial Estate, 9 LBS Marg, Kurla (W), Mumbai 400070
Tel: (022) 4215 5441 / 93204 22266

Hyderabad: H.NO-3-41-66, Sri Laxmi Nilayam, South Block, West Marredpally, Secunderabad 500026
Tel: 93468 32266

Nagpur: 93738 38883, **Ahmedabad:** 93284 22333, **Lucknow:** 93361 20941

E-mail: sales@nidekindia.com, Website: www.nidekindia.com

A product of
Nidek Medical Products Inc.
Birmingham
Alabama 35217, USA