



REVEALA

3 Channel Digital Holter

A Holter in rhythm with your practice

Key Features

- ✓ Available in 3 & 12 channel models.
- ✓ Detects 12 types of arrhythmia.
- ✓ Heart Rate Variability (HRV) analysis.

REVEALA 3 Channel Digital Holter



Unit Specifications

Instrument Type	: 3 /12 channel solid state Holter.
Acquisition	: Simultaneous acquisition of leads.
Recording time	: 24 Hrs. max, continuous non-stop.
Indications	: Battery Low (audio visual) Leads off (audio /visual) Memory full (visual) Recording status (visual) Event mark (audio)
Host interface	: USB interface
Storage media type	: in-built SD card
Power requirement	: Powered by (1.5v) alkaline batteries 2 AA size (3ch.), 2 AA size (12ch.)

Analysis software specifications


Configurable disclosure time of 10, 20, 30, 60 minutes.

Selection of any of the channels for disclosure.

Review parameters	: Min HR. Longest pause Longest SVTach. Max delta ST in mm with channel and time. Heart rate variability (HRV) analysis.	: Max HR. Longest Vtach. Longest ST episode.
Trends	: HR V tach SVE tach Trigeminy	: VE SVE Bigeminy Bradycardia, Tachycardia
Arrhythmia	: Auto detection Pause Trigeminy VE pair V-Tach SV pair Bradycardia	: Bigeminy Isolated VE R on T SVE SV Tach Tachycardia
Print Reports	: Arrhythmia disclosure & summary Full disclosure-1 min to 60 min (user selectable) Template summary ST Analysis report Summary report. Heart rate variability report	

Minimum PC configuration

Operating system	: Win7 / Windows XP / Windows 2000 P4, 2 Ghz, 20 Gb HDD 256 MB RAM, USB port Super VGA colour monitor
------------------	--

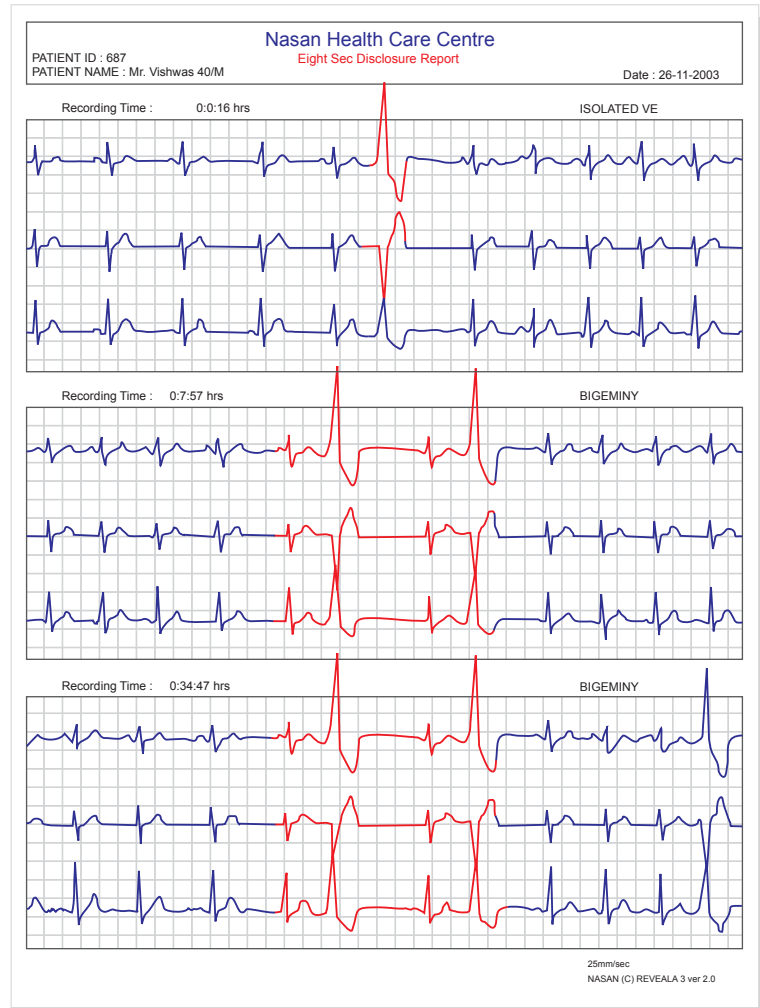

DUE TO OUR CONTINUOUS PRODUCT DEVELOPMENT PROGRAM, THE SPECIFICATIONS MAY CHANGE WITHOUT PRIOR NOTICE.

NASAN MEDICAL ELECTRONICS PVT. LTD.

Plot No. 147, "Rangali" Society, First Floor, Opp. Karve Statue, Pune - 411 038.

Phone : +91 82370 32389 Fax : (020) 25393540

E-mail : sales@nasanmedical.com Website : www.nasanmedical.com